


New format

**PRODUCT INFORMATION CONCERNING
PLATELET CALIBRATOR KIT (REF 7011/00457) – NEW FORMAT**

Dear "PLATELET Calibrator" kit users,

We would like to inform you that the "PLATELET Calibrator" kit features now a new calibrator reagent.

In order to broaden the use of the "PLATELET Calibrator" kit to newer cytometers bearing different optical properties, we have decided to switch from a calibrator made of 4 populations of beads to a calibrator consisting of 3 populations of beads. This modification does not change the performances of the kit.

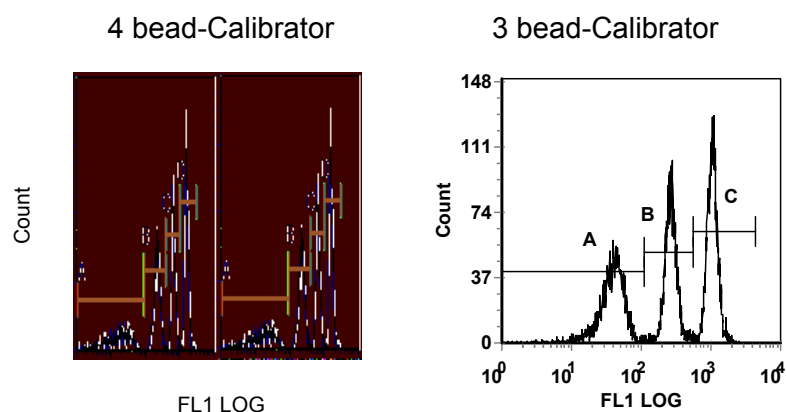

Indeed we have noticed that the lowest calibration bead could generate slightly different quantitation values between cytometers from older and newer generation.

This new calibrator will lead to homogenous quantitative results between cytometers from different generation and allow to switch the analysis protocol to a newest instrument without any difficulty.

At the first analysis of this new calibrator, check the PMT FL1 in order to get, on the FL1 Log vs Count histogram, the third calibration bead at the beginning of the 4th decade and adjust the cursors on the 3 peaks of fluorescence.

Once this optimisation is done, you will analyse your samples as usual.

Yours sincerely,

The Biocytex Marketing Department

Distributed by:



MILAN ANALYTICA AG
YOUR BIOTECHNOLOGY PARTNER

Baslerstrasse 15
4310 Rheinfelden
Switzerland
www.milananalytica.ch
info@milananalytica.ch

Tel.: +41 61 845 99 88
Fax: +41 61 845 99 87